



CRYO-NUMBING DEVICE

CRYO TECHNOLOGY

NO CHEMICALS

NUMBS THE NEURONS AT THE SITE OF APPLICATION

ONLY TAKES UP TO 10 SECONDS TO NUMB

NUMBS AT -9°C

INJECT IMMEDIATELY AFTER

PRIORITISES PATIENT COMFORT

COMPACT IN SIZE

EASY TO READ LIGHT INDICATOR

1500+ USES PER DEVICE



EXOCOOL

Just freeze, apply and inject! EXOCOOL is a non-invasive cryo-numbing device that is designed to make injections painless. The cryo-numbing mechanism numbs the area of application to -9°C in just 10 seconds. Made from medical-grade ABS plastic and aluminium alloy, each EXOCOOL device gives up to 1,500 uses. Excellent for cosmetic procedures and is safer and faster than numbing cream and ice packs.

APPLICATIONS

Use EXOCOOL to numb the area for:

- Lip Filler
- Botox
- Facial Filler
- Filler Dissolving
- Scalp regeneration
- Multi-needle injections
- PRP
- Extractions
- Non-surgical soft blepharoplasty
- Tattoo removal
- Threads
- Eyebrow Waxing, Plucking & Threading



User Friendly

Hand-held, convenient, and easy to use.



Multiple Applications

One solution for all kinds of needle pain.



Targeted Application

Works effectively on the site of application.



No Reaction

Non-invasive i.e no chemicals used



Quick Action


Numbs the area within 6-10 seconds.



Enhances Comfort


Makes every kind of injection/needle painless

EASY TO READ LIGHT INDICATOR

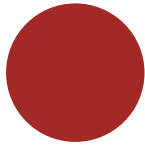

Red Light

The Device isn't cold enough.
Leave in the freezer for a longer period of time

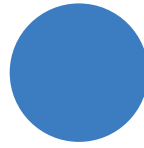

Blue Light

The Device is too cold. Leave outside the freezer for a few minutes



Green Light

The Device is at optimum temperature and is ready to use.

SAFE, CONVENIENT AND EFFICIENT

EXOCOOL	TRADITIONAL NUMBING CREAM
Chemical Free	Has Chemicals
No Mess Application	Messy application
Only takes 4-6 seconds to numb sensitive skin & Only takes 8-10 seconds to numb normal skin	Can take up to 40 minutes to numb
Can be used around the eyes by putting a rubber glove around the metal tip instead of using the Gel	Can't use around the eyes due to reaction
More cost efficient with 1500 uses per device	Requires constant repurchasing with risk of expiration and need for prescription / Dr approval
No Adverse reactions	Numbing cream has risk of adverse reactions

EXOCOOL	GENERIC ICE PACK
Numbs the area at -9°C	Numbs at only 0°C
Easy and more comfortable to hold as the metal tip is the only area of the device that is cold	Uncomfortable for the patient and practitioner to hold
Easy to read indicator light to determine if the device is at the optimum numbing temperature	No clear indicator
Has targeted application with pinpoint accuracy	No pinpoint accuracy, where it numbs all surrounding areas
Localised Application	Unsuitable for small areas

3 SIMPLE STEPS



Step 1: Freeze

Freeze the device by keeping it in the deep freezer at a temperature below -12°C .



Step 2: Apply

Check if the light is blinking green. If it is, apply the metal tip on the required area for 4-10 seconds*



Step 3: Inject

Administer the injection to the area of numbing for a painless injection experience

A Revolutionary Pain-numbing Device

With Exocool, say goodbye to the fear and discomfort of the sharp prick of injection needles, and other cosmetic procedures and treatments.



Product Size	160mm x 44mm x 36mm
Product Weight	120g
Operating Temperature	-20°C to 50°C
Duration of Exocool Application	4-6 Sec for Sensitive Skin 8-10 Sec for Normal Skin
Light Indication (Usage Indication)	Blinking Green, Blinking Red, Blinking Blue
Battery Capacity	2500+ Checks
Product Life	1500 Uses, Under Ideal Conditions
Product Material	Virgin Medical Grade ABS-plastic and Medical Grade Aluminium Alloy
Cooling Reduction Rate (with Cap)	0.5°C per Minute at Room Temperature (27°C)
Cooling Reduction Rate (without Cap)	0.7°C per Minute at Room Temperature (27°C)
Advisable Cooling Temperature	Below -12°C
Minimum Cool Time (After Red Light)	30-45 Minutes, Below -12°C